



**HUMAN CARE**  
Your Life. Your Way.



**Convertible Chairs**

# Introduction

We produce and provide lifting solutions for people with special needs. Our products are made to serve as a natural part of life to all our users.

Human Care has a proud history in patient care. Founded in 1993, our twenty plus years of experience is expressed through a broad selection of premium products offering high customer value and satisfying all customer needs.

Our main product areas are: Mobility Aids, Lifting Solutions, Healthcare Beds, Convertible Chairs, and Bathroom Safety Products for people with reduced mobility.

With headquarters in Stockholm, Sweden, and subsidiaries in the USA, Canada the Netherlands and Australia, we collaborate with distributors to provide solutions in more than twenty countries.

Our promise is safety, quality and usability. In all of our products, we exercise only the highest quality to make sure that every detail, every aspect, every function will deliver that promise made to our users.

**Welcome to Human Care.**



# Content

	Page
H-250	7
I-400	8
I-700	9
Patient Transfer System	10

# Effective and safe early mobilization with the Human Care Convertible Chair

The Human Care Convertible Chairs offer unlimited number of reclined and tilt positioning. The ability to frequently reposition a user significantly increases the users tolerance for long-term seating. A sturdy frame can be equipped with a patient transfer system to enable a single caregiver to safely transfer a user to and from bed.

## Clinical benefits for the patient with the Human Care Convertible Chair

### Faster recovery

The ability to get the patient out of bed faster can reduce the length of stay in both higher acuity departments and the overall length of stay in the hospital.

### Safe and effective transport

Transferring a patient in a convertible chair will facilitate transfer of patients to allow transport to physical therapy and/or other ancillary services.

### Patient experience

Not every patient can be lifted. Some patients have pressure ulcers, pain or are afraid to be lifted. The convertible chair is a non-lifting alternative that will create the best transfer possible for these patient groups.

## Benefits for the caregiver with the Human Care Convertible Chair

### Makes your job easier

The convertible chair makes your job easier and helps enhance the quality of care provided while reducing recovery time. It has been proven to deliver confidence to the patient and their loved ones to facilitate early mobilization and rehabilitation.

### Efficient and easy to use

Once transferred to the convertible chair a patient can easily be repositioned by simply using all the unlimited number of tilting options.

### Injury prevention/

### Ability to drives SPH initiatives

Facilities utilize the convertible chair to protect the clinicians from injury and to help drive SPH initiatives.







” Our mission is to improve the quality of care and freedom of movement for people with special needs

# Investment in an ICU early rehabilitation program can generate net financial savings for hospitals

## Financial Benefits with the Human Care Convertible Chair

### **Reduced length of stay**

Studies have shown that early mobilization leads to reduce average length of stay in the hospital for a give episode of care.

### **Reduce cost of patients in high acuity departments**

Studies have shown that early mobilization leads to lower costs per patient in departments that require getting a patient to rehab and into step down units.

### **Helps to reduce injuries to caregivers**


Reiterate cost to hire a new clinician/nurse, ST and LT costs to replace a clinician.

### **Supports/drives SPH initiatives**

Complimentary to investments to existing initiatives/solutions. Chair can also help facility to initiate SPH in highest cost/risk departments.

### **Efficiency**

Ability to utilize one solution for patient's room.

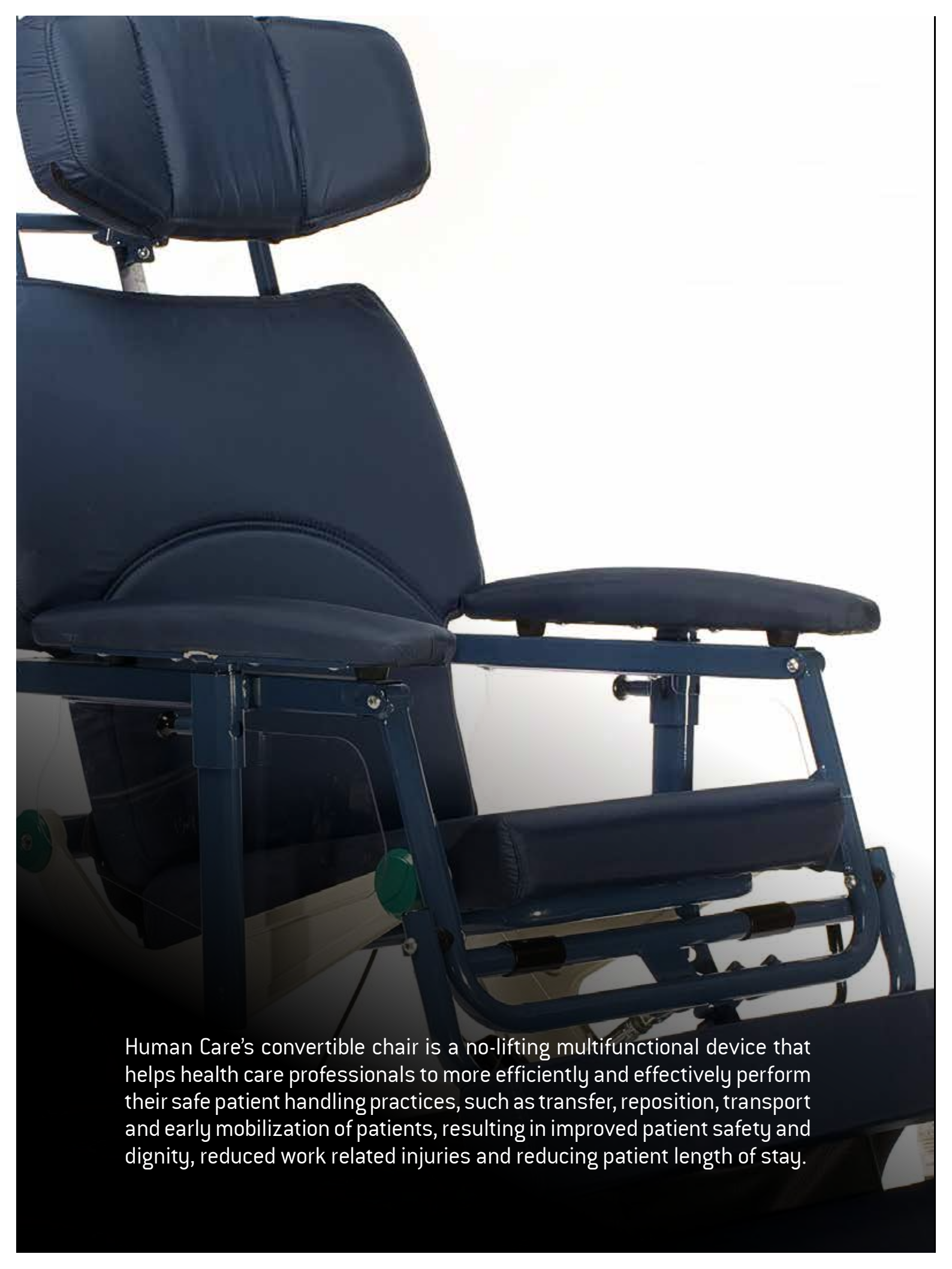

Human Care's convertible chair is a no-lifting multifunctional device that helps health care professionals to more efficiently and effectively perform their safe patient handling practices, such as transfer, reposition, transport and early mobilization of patients, resulting in improved patient safety and dignity, reduced work related injuries and reducing patient length of stay.

## H-250

Convertible chair for moving, repositioning and transferring patients

- Provides safe transfer to and from bed
- Limitless number of positions
- Ideal for home use



The H-250 Convertible Chair is a unique device that easily converts from a stretcher to an upright mobile chair in seconds. Its Tilt-in-Space capability allows for increased blood circulation helping to prevent skin breakdown and other complications of immobility.

With its **113 kg/250 lbs capacity**, the H-250 provides a no-lift solution for moving, repositioning and transferring patients in the safest manner possible.

The H-250 chair together with a Patient Transfer System (PTS) offers the ability to move between a bed and a chair with ease and enables a single caregiver to easily transfer a patient. Read more about our PTS on page 10.



Description	Part no.	Seat width	Overall width	Headrest adjustable height	Maximum user weight	Accommodates varying patient heights
H-250 Chair	B1002-BL	431 mm / 17"	610 mm / 24"	178 mm / 7"	113 kg / 250 lbs	1270 - 1980 mm / 5' 0" - 6' 6"

## I-400

### Convertible chair for moving, repositioning and transferring patients

- Unlimited repositioning capabilities
- Easy to turn and maneuver
- Increases the users tolerance for long-term seating



The I-400 chair is a convertible chair for safe and comfortable moving and transferring of patients. It has unlimited number of reclined and tilt positions. The ability to frequently reposition a user significantly reduces the chance of pressure ulcers and increases the users tolerance for long-term seating.

The I-400 chair together with a Patient Transfer System (PTS) offers the ability to move between a bed and a chair with ease and enables a single caregiver to easily transfer a patient. Read more about our Patient Transfer System on page 10.



Central locking/steering casters safe for patient and caregiver



Can be equipped with a patient transfer system. See page 10.



Product	Part no.	Seat width	Overall width	Overall length	Maximum user weight	Accommodates varying patient heights
I-400 Chair	B1041-TS-GEL	546 mm / 21.5"	749 mm / 29.5"	1880 mm / 74"	181 kg / 400 lbs	Manual Tilt and Recline



# I-700

Convertible chair for moving, repositioning and transferring patients

- Tilt-in-space function
- Helps to reduce risk of injuries
- Fully adjustable head support



The I-700 chair offers an unlimited number of reclined and tilt positions. To frequently reposition a user significantly reduces the chance of pressure ulcers and increases the users tolerance for long-term seating.

The I-700 chair can be equipped with a patient transfer system to enable a single caregiver to perform safe lateral transfer to and from bed. Read more about our Patient Transfer System on page 10.



Fully adjustable headrest with support



Can be equipped with a patient transfer system. See page 10.



Product	Part no.	Seat width	Overall width	Overall length	Maximum user weight	Accommodates varying patient heights
I-700 Chair	B1071-TS-GEL	686 mm / 27"	889 mm / 35"	1880 mm / 74"	318 kg / 700 lbs	Manual Tilt and Recline

# PTS

## Patient Transfer System

- Allows a single caregiver to safely transfer a patient
- Provide the safest transfer possible
- Create a more positive patient experience



Moving a patient can be challenging and unsafe, especially when dealing with heavier patients. All our chairs can be equipped with our unique, easy to adjust Patient Transfer System (PTS).

**Our Patient Transfer System works with** most beds and stretchers and allows a single caregiver to safely and easily transfer a dependent patient with greater frequency and ease without the indignity, discomfort and medical limitations of traditional handling methods.



Product	Included in kit	Part no.
Patient transfer system kit for H-250	2 chair brackets, 5 transfer straps, 2 bed posts, 1 mattress cover (non-slip), 1 Transfer bar, 5 sheet clips, 5 plastic plugs and 2 bed clamps.	b2004e
Patient Transfer System kit for I-400 (hospital)	2 bed brackets, 1 transfer bar, 5 transfer straps and 2 transfer sheets.	b2011a
Patient Transfer System kit for I-400 (Home care)	2 bed posts, 1 transfer bar, 5 transfer straps, 1 mattress cover (non-slip), 5 plastic plugs, 5 sheet clips, 2 bed clamps.	b2001c
Patient Transfer System kit for I-700	2 bed brackets, 1 transfer bar, 5 transfer straps and 2 transfer sheets.	b2013a

# Patient Transfer System – Easy transfer performed by one single caregiver

The duties today required of caregivers have become increasingly dangerous due to the constant demand of moving and repositioning patients. This is especially true working with bariatric patients.

Through our PTS system (Patient Transfer System) unnecessary lifting and risk of injury to caregivers is minimized. Our Patient Transfer System (PTS) together with our convertible chairs offer a multitude of positioning capabilities within a single product.


Designed to allow a single caregiver to easily and safely transport and reposition a patient, our chairs can perform the functions of stretcher, mobile chair, and custom seating.

## The Benefits of the Patient Transfer System

The work required of caregivers has become increasingly dangerous due to the constant demand to move and transfer dependent patients. So, we're determined to:

- Make this job easier
- Reduce preventable incidents
- Provide the safest transfer possible
- Help enhance the quality of care provided
- Create a more positive patient experience


**1** Place the chair by bed in stretcher position.

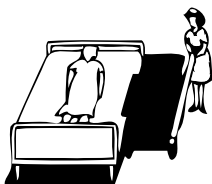


**2** Attach bed sheet to transfer bar.



**3** Transfer patient by turning handle.

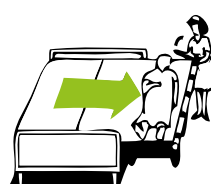

**4** Incline chair and you're ready to go!

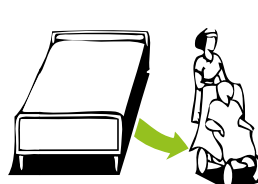



Distributor:



**HUMAN CARE**  
Your Life. Your Way.

**AUSTRALIA**

Unit 4, 49 Bertie Street  
Port Melbourne, VIC 3207  
Phone : +61 1300 132 736  
Fax : + 61 [03] 9645 4954

[info.au@humancaregroup.com](mailto:info.au@humancaregroup.com)

**CANADA**

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: 613.723.6734  
Fax: 613.723.1058

[info.ca@humancaregroup.com](mailto:info.ca@humancaregroup.com)

**SWEDEN**

Årstaängsvägen 21C  
117 43 Stockholm  
Phone: +46 8 665 35 00  
Fax: +46 8 665 35 10

[info.se@humancaregroup.com](mailto:info.se@humancaregroup.com)

**NEDERLAND**

Elspeterweg 124  
8076 PA Vierhouten  
Phone: 0577-412171  
Fax: 0577-412170

[info.nl@humancaregroup.com](mailto:info.nl@humancaregroup.com)

**UNITED STATES**

4210 S. Industrial Drive, Suite 160  
Austin, TX 78744  
Phone: 512.476.7199  
Fax: 512.476.7190

[info.us@humancaregroup.com](mailto:info.us@humancaregroup.com)