

The Mobi Pro Flexi is an ergonomic and compact mobile assistive mobility device ideal for people with low leg strength and some trunk balance who experience difficulty walking. The Mobi Pro Flexi is an excellent safe patient handling solution for early mobilization. It offers support while standing and is an efficient device for patient transfers, improving safety and efficiency.

#### Rotating Seat

The Mobi Pro Flexi is an efficient and safe tool for the healthcare provider to move a patient, improving safety and caregiver time management. It features rotating padded seats and knee support for easy application. The feet are placed on a solid, firm footplate and then the patient pulls himself/herself up to a standing position using the integrated solid bracket. By simply rotating the seat cushion, the Mobi Pro becomes a comfortable seat that the patient can rest and be moved on.

#### Versatile Design

Using the Mobi Pro Flexi not only encourages mobilization and activity but gives the caregiver fast, stable and safe transport of a

patient to the bed, chair or the toilet and back if needed. The Mobi Pro Flexi is equipped with a mechanical mechanism to spread the legs, ensuring it will work with most seating (like a toilet or wheelchair). The vast adjustment range is easily set to three settings with a lever.

#### Quality

Our products are designed and manufactured to meet all the requirements of relevant standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone: +46 8 665 35 00  
Fax: +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street  
Port Melbourne, VIC 3207  
Phone: +61 1300 132 736  
Fax: + 61 (03) 9645 4954  
info.au@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com

#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories



**Foot Plate**  
Included with device.  
Durable ABS plastic.  
Replacement Part no.  
1590035



**Casters**  
Double wheel 7.5cm/3".  
Front Part no. 1690060  
Rear with Brake Part no.  
1690061



**Ergonomic Grips**  
Made of rubber. Sold in  
pairs. Replacement Part  
no. 1510054

## Lifting Situations



**Use to/from a bed**  
Use to raise/lower a  
patient from/to a bed.



**Use to/from a chair**  
Use to raise/lower a  
patient from/to a chair.



**Use to/from the toilet**  
Use to raise/lower a  
patient from/to a toilet.




Tested to all relevant  
standards, including ISO  
14971, EN12182, MDD  
93/42/EEG.

## Features



**Easy Movement**  
The Mobi Pro Flexi is  
easily moved by a single  
caregiver. A small turning  
radius makes it ideal for  
smaller spaces.

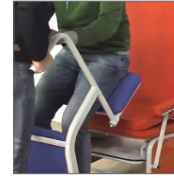

**Adjustable Base**  
Manual adjustment of  
legs allows it to fit under  
most seating.



**Padded Rotating Seats**  
Allows device to be  
used for seating as  
an alternative to a  
wheelchair.



**Padded Knee Pad**  
Soft, cushioned knee pad  
to add patient comfort  
and reduce device  
pressure at the knees  
while standing.

## Technical Specifications

Mobi Pro Flexi			
Part No.	92500	Adjustable Leg (Min - Max)	43cm - 70cm / 17" - 27.5"
Width	58.5cm / 23"	Handle Height (Inclined at Angle)	92cm - 100cm / 36.2" - 39.4"
Length	97cm / 38.2"	Weight	19.6kg / 43.2lb
Height	104cm / 41"	Weight Capacity	150kg / 330.6lb
Wheel Diameter	7.5cm / 3"		

This product is manufactured  
in accordance with all  
relevant ISO standards.



Learn More  
About This Product



### Follow Us on Social Media



facebook.com/humancarecda



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fyjNH3zTjRxOP1IDf2DaQ